



Designed to Last





Design

The VIDAS® range is greatly
Even with years of experie
So, after listening to your cu
we've taken this already in



Unique values with a stylish-design puts VIDAS

Improved Design for Increased User Comfort

From your requirements we have advanced the system's design to further :

- **Streamline your workflow** by optimising test loading
- **Increase user-friendliness** for data entry and reporting patient results.

Timeless VIDAS Values

- Single-dose, ready-to-use reagents (1 patient = 1 test)
- All inclusive reagent kits
- System robustness (MTBF* > 452 days)
- Calibration every 14 days
- Automated barcode identification
- Just Load and Go

Flexibility, where and when you need it :

- Five independent 6-test sections
- Processing of different parameters simultaneously
- Large routine menu (over 80 parameters in Immunochemistry and Infectious Diseases)
- Packaging suited to different test volumes (30 or 60 tests)

Quality

- No cross-contamination
- Reference tests : HIV DUO Ultra, D-Dimer Exclusion™, Toxo Avidity...
- CVs highlighted as being excellent by numerous external quality controls

* Mean Time Between Failure

VIDAS Size:


58 cm (23 inches) (height)
105 cm (41.5 inches) (width)
68 cm (21.75 inches) (depth)

VIDAS weight:

65 kg (143 lbs)

Electrical requirements:

Voltage: 100-240 VAC
Electrical Consumption : 3-1.2 A
Frequency: 50/60Hz
Power: 280 VA maximum
Heat emission: approx.921 BTU/hr

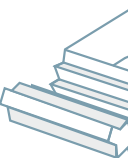

REAGENT

d to Last

appreciated worldwide for its simplicity, flexibility and accuracy.
ce, our continued success relies on you.
urrent needs to further optimise workflow
novative design one step further.

AS in a class of its own.



Simplicity makes the difference


The VIDAS test system uses a single-dose, ready-to-use reagent.

- 1• The **SPR**[®] acts as a Solid Phase Receptacle for the reaction and is coated with antigens or antibodies.
- 2• The **Strip** contains all the reagents required for the reaction.

The SPR acts as a pipetting device.

At each stage of the reaction, it aspirates the reagents in and out. This original concept prevents any inter-reagent or inter-sample contamination. In addition, the absence of tubing, syringes and needles reduces system maintenance to a minimum.

The VIDAS system uses the ELFA assay principle, combining the ELISA test method with a final blue fluorescent reading. This technology ensures excellent result sensitivity and specificity.

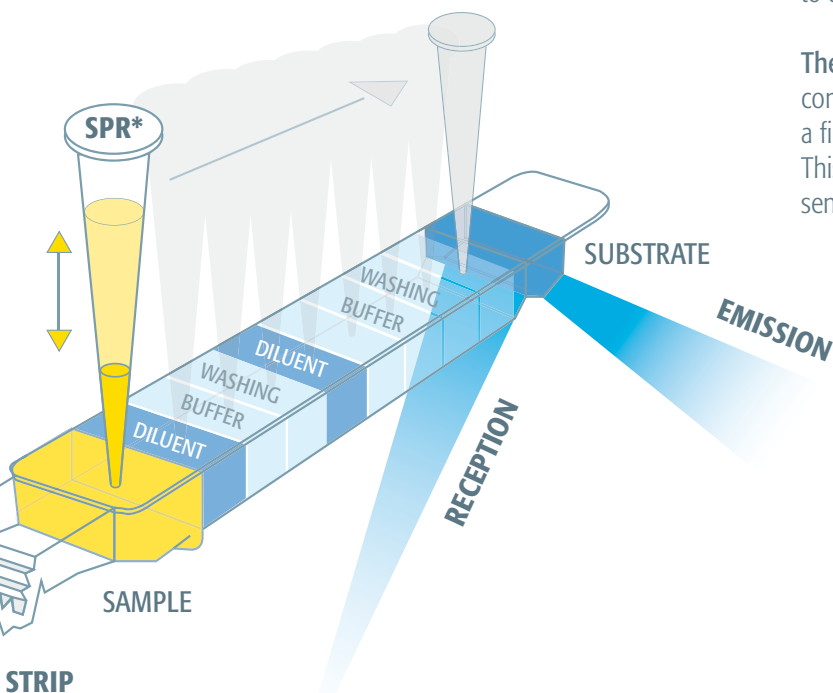

* SPR: Solid Phase Receptacle

VIDAS, Designed to last

Ease-of-use: Load and Go

Robustness: MTBF > 700 days

High Quality tests for reliable diagnosis

Cost effective solution

- **Single dose**
- **All inclusive kits**
- **Convenient packaging size**

VIDAS and mini VIDAS : adapted offer to your organisation.

VIDAS, a Customized Solution for your Daily Workload

Routine Testing

- Up to 10 parameters in one run
- Comprehensive panels (Thyroid, fertility, cardiac markers, tumor markers, HIV, hepatitis, ToRC...)

Complementary and Specialized Testing

- Confirmatory testing (HIV DUO, Toxo Avidity, DD Exclusion, Prolactin...)
- Specialized testing: Lyme, Measles, Mumps, H. pylori, CD A/B, Chlamydia, Allergy...)
- Complementary Testing: large menu (> 90 parameters)

Emergency Testing

- 24/7 available
- Rapid results in compliance with Healthcare guidelines
- Emergency panel: TNIU, CK-MB, Myoglobin, NT-proBNP, DD Exclusion, PCT, Digoxin, hCG, p24 Ag...



bioMérieux S.A.

69280 Marcy l'Etoile

France

Tel.: 33 (0)4 78 87 20 00

Fax: 33 (0)4 78 87 20 90

www.biomerieux.com

www.biomerieux-diagnostics.com

