


Maximize Your Efficiency With Myla® Microbiology Middleware Solution



Optional advanced reporting and connectivity with Myla®

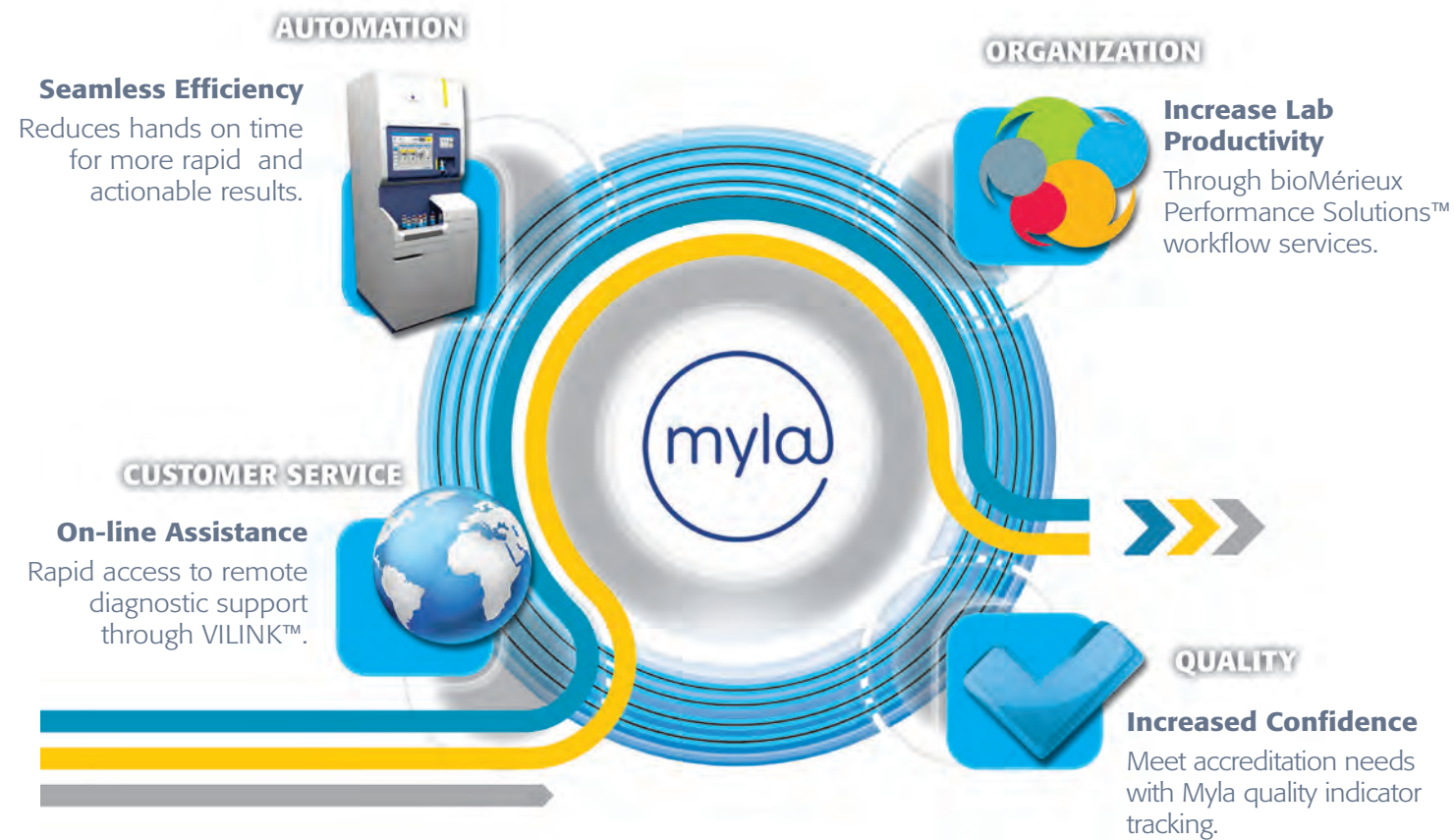

Connect multiple instruments and gain a comprehensive view of microbiology workflow with Myla. Myla provides quick access to critical information needed to assess the overall quality of the blood culture testing process:

- Contamination tracking
- Blood volume quality indicator reporting
- Time to detection for blood culture isolates
- Positivity and negativity rate analysis

04-14 / 9307402/002/CB/A / Document and/or picture not legally binding. Modifications by bioMérieux can be made without prior notice. / BIOMERIEUX, the blue logo, BIOMERIEUX PERFORMANCE SOLUTIONS, FAN, MYLA, BacT/ALERT, VILINK and VIRTUO are used, pending and/or registered trademarks belonging to bioMérieux, or one of its subsidiaries, or one of its companies. Any other name or trademark is the property of its respective owner / bioMérieux S.A. RCS Lyon 673 620 399 / Photos : 3Dvision Nyalls / THERA Conseil / Printed in France.

BioMérieux S.A.
69280 Marcy l'Etoile
France
Tel. : 33 (0) 4 78 87 20 00
Fax : 33 (0) 4 78 87 20 90
www.biomerieux.com
www.biomerieux-diagnostics.com

Manufactured by :
bioMérieux Inc.
100 Rodolphe Street,
Durham, NC 27712
USA
bioMérieux Inc. is a wholly-owned subsidiary of bioMérieux sa.

THE NEXT DIMENSION IN BLOOD CULTURE DETECTION

THE NEXT DIMENSION IN BLOOD CULTURE DETECTION

VIRTUO™, refining and redesigning your workflow for maximum laboratory efficiency.

MOTION ACTIVATED LOADING

Motion activated loading bay with optical sensors.

Detects user intervention.
Conveyor and carousel feed BacT/ALERT bottles into instrument.

SMART SCANNING TECHNOLOGY

Carousel moves BacT/ALERT bottles into **smart scanning station**.
Scanning station rotates and **images entire BacT/ALERT bottle**.

SAMPLE LEVEL SENSING

Directly senses fill level **on every BacT/ALERT bottle** in real time.

PRECISION LOADING

Robotic arm grips BacT/ALERT bottles and places firmly into cells.
Automatically performs periodic **QC on each cell**.

SYSTEM PERFORMANCE

New growth detection algorithms for overall **faster time to detection**.

AUTOMATED UNLOADING

Final **negative BacT/ALERT bottles removed automatically** or on demand.
Shatter resistant BacT/ALERT bottles for worry free unloading.



Seamlessly Efficient blood culture technology in one automated system to create a new dimension in blood culture detection.

VIRTUO™ seamlessly integrates your blood culture workflow to provide efficient and actionable results.

- Motion activated "Set and Forget" loading so any skill level can load BacT/ALERT® bottles.
- Intuitive touchscreen design for enhanced ease of use.
- Automated unloading with removable waste container.
- VIRTUO helps you redesign your workflow by reducing hands on time and increasing efficiency.



THE NEXT
DIMENSION
IN BLOOD
CULTURE
DETECTION